

# EFFICIENT HYDROPHILIC INTERACTION CHROMATOGRAPHY [HILIC] METHOD DEVELOPMENT STRATEGY




## ACQUITY UPLC Columns and XBridge HPLC Columns

### Ethylene Bridged Hybrid [BEH] Phases

#### BEH HILIC

- Unbonded Ethylene Bridged Hybrid (BEH) substrate
- Particle Size: 1.7 µm [UPLC]; 2.5, 3.5 and 5 µm [HPLC]
- Pore Size: 130Å; Surface Area: 185 m<sup>2</sup>/g; Carbon Load: NA
- pH Range: 1-9; Low pH Temp. Limit: 45 °C; High pH Temp Limit: 45 °C
- Recommended Usage: Excellent for the retention of very polar, basic, water soluble analytes. Specifically designed and tested for HILIC separations using high organic mobile phases.

#### BEH Amide

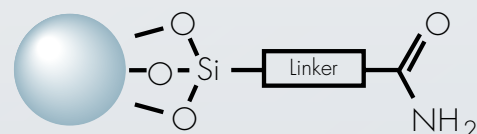

- Trifunctional amide, bonded to Ethylene Bridged Hybrid (BEH) substrate
- Particle Size: 1.7 µm [UPLC]; 3.5 µm [HPLC]
- Pore Size: 130Å; Surface Area: 185 m<sup>2</sup>/g; Carbon Load: 12%
- pH Range: 2-11; Low pH Temp. Limit: 90 °C; High pH Temp Limit: 90 °C
- Recommended Usage: Rugged HILIC stationary phase designed to separate polar acidic, neutral and basic compounds. Particularly well suited for carbohydrate [saccharide] analysis using high concentrations of organic modifier, elevated temperature and high pH.



## Atlantis HPLC Columns

#### Atlantis HILIC Silica

- Unbonded high purity silica substrate
- Particle Size: 3, 5 and 10 µm [HPLC]
- Pore Size: 100Å; Surface Area: 330 m<sup>2</sup>/g; Carbon Load: NA
- pH Range: 1-5; Low pH Temp. Limit: 45 °C; High pH Temp Limit: 45 °C
- Recommended Usage: Excellent for the retention of very polar, basic, water soluble analytes. Specifically designed and tested for HILIC separations using high organic mobile phases.

## ESSENTIAL INFORMATION FOR METHOD DEVELOPMENT

### Sample solubility

- A mixture of 75% acetonitrile with 25% methanol is a good general purpose diluent for the best combination of sample solubility and peak shape. 0.2% formic acid or ammonium hydroxide can be added to further encourage analyte solubility.
- The best chromatographic performance will be observed when the sample is dissolved in the initial mobile phase composition.
- Peak distortion may occur when the sample is dissolved in a solvent of very different polarity and viscosity to that of the mobile phase.
- The impact of varying the injection volume should be investigated.

### Chemical properties: functional groups

- Analytes will often exhibit the most retention when they are ionized [i.e., bases at low to mid pH, acids at mid to high pH], which is the opposite of reversed phase.
- For sample mixtures containing different types of analytes, use an intermediate pH ~ 5-7 [10 mM ammonium acetate with 0.02% acetic acid].
- Knowledge of the pKa value[s] of the analytes is helpful when developing a method.

### Mobile phase preparation

- For the best gradient performance and reproducibility, a constant ionic strength of buffer [10 mM] or additive [0.2%] is recommended.


### Column equilibration

- Condition new columns with 50 column volumes of 50:50 acetonitrile:water with 10 mM buffer or 0.2% additive solution. Store the column in these conditions.
- Equilibrate the column with 20 column volumes of initial mobile phase conditions before routine use [approximately 2 blank gradients].
- For gradient analysis, re-equilibrate with 5-8 column volumes between injections.


### Needle wash solvents

- The needle wash solvent should closely match the initial mobile phase conditions [90-95% acetonitrile] to preserve retention and peak shape.


## SYSTEMATIC SCREENING METHOD DEVELOPMENT PROTOCOL

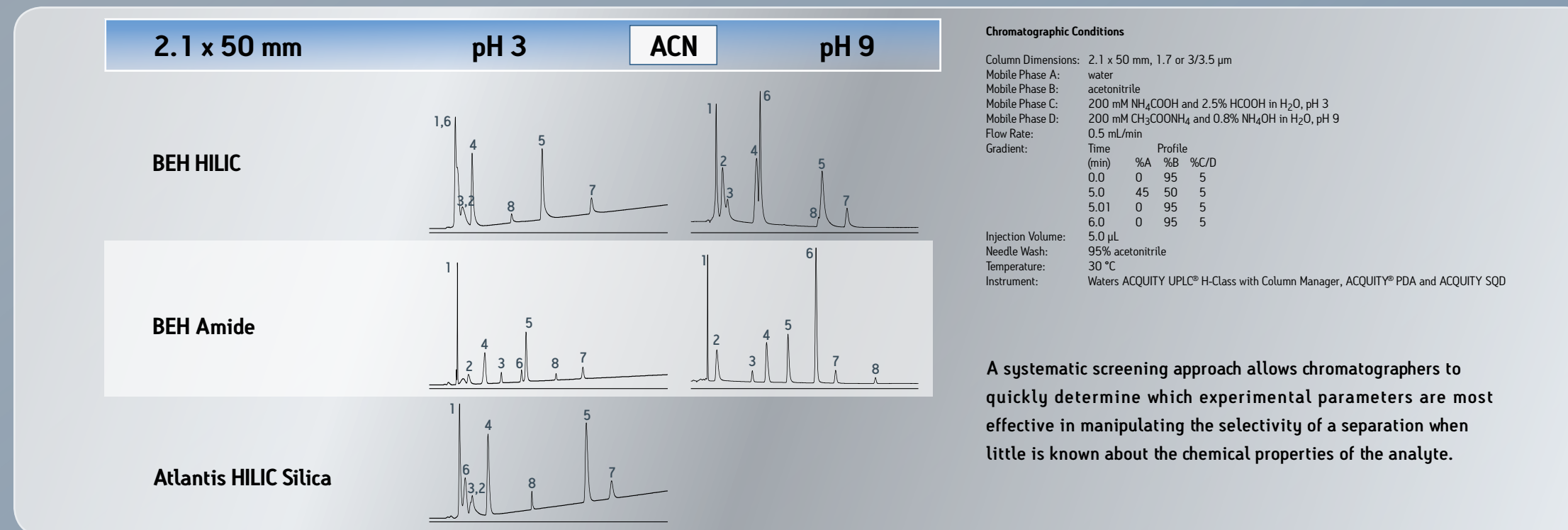


## THE IMPORTANCE OF MOBILE PHASE pH

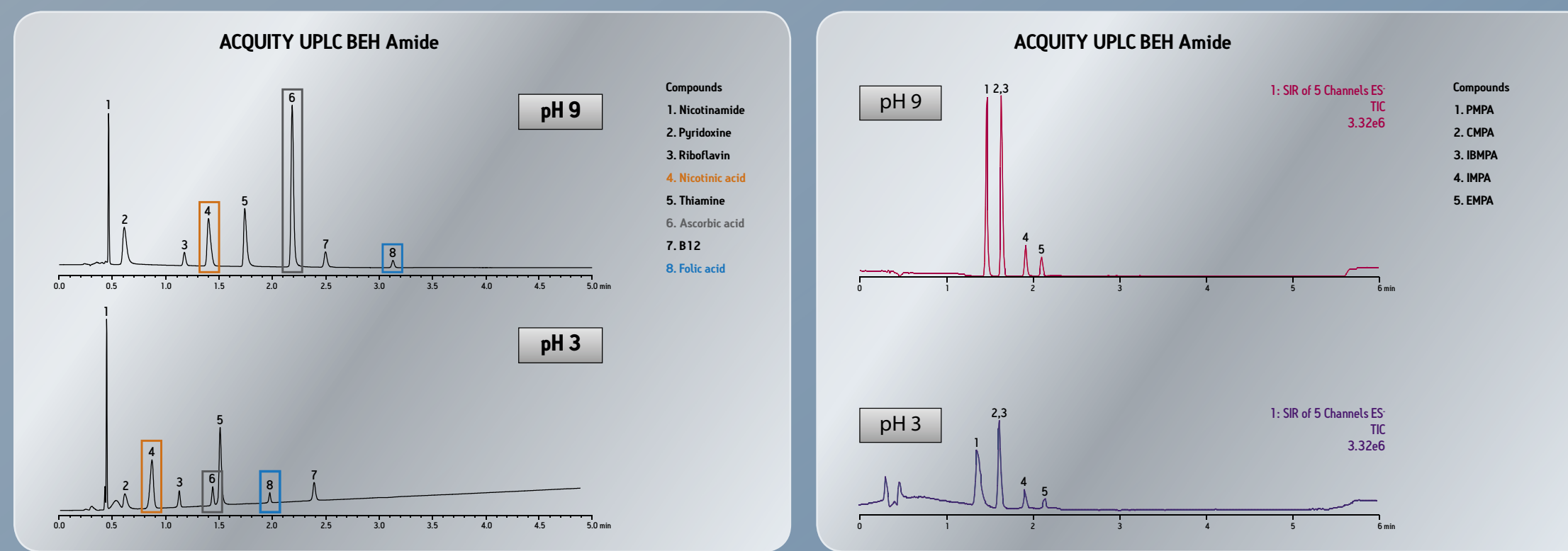

## ORGANIC MODIFIER SELECTION

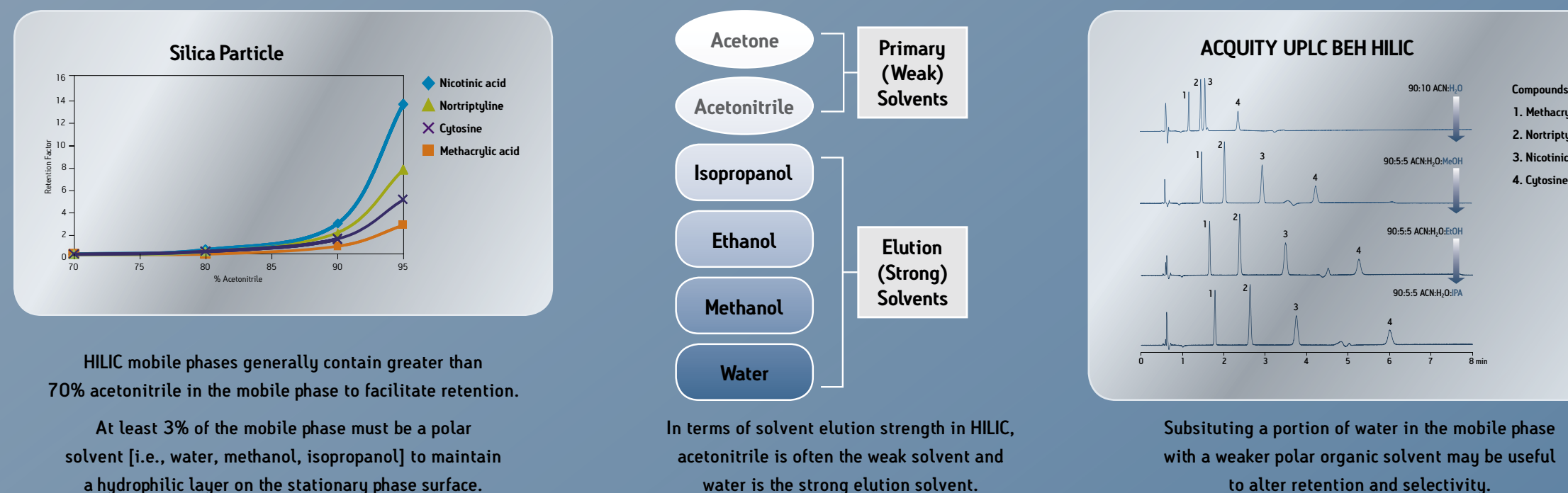


## METHOD SELECTION AND OPTIMIZATION

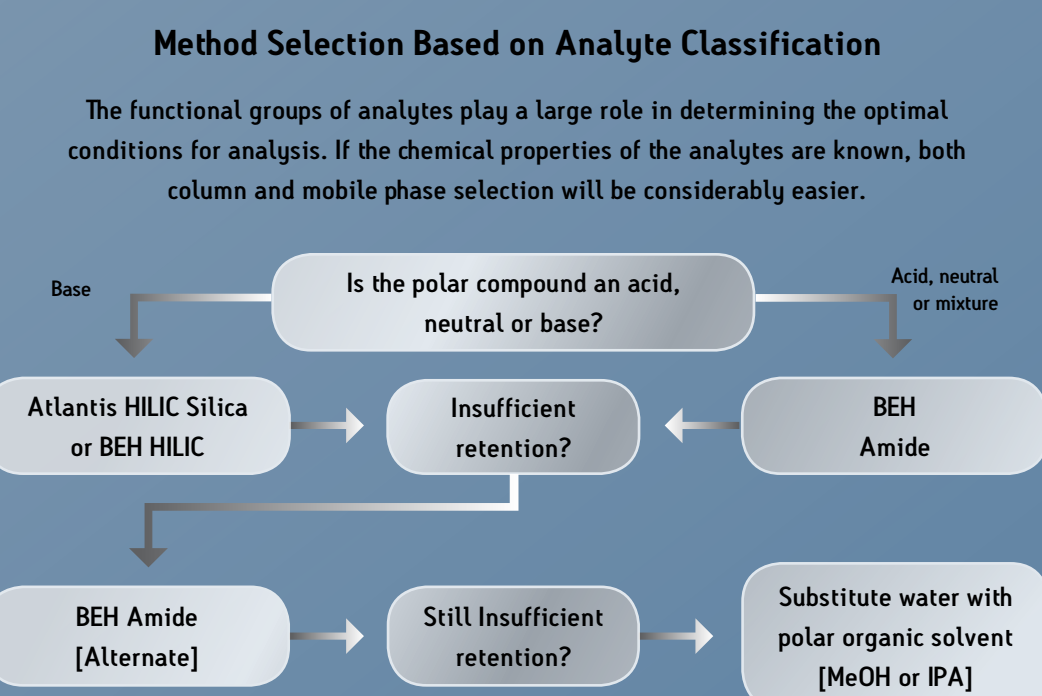

**Key**  
 Atlantis HILIC Silica = low pH (pH 3), 10 mM ammonium formate with 0.2% formic acid in 95:5 ACN:H<sub>2</sub>O  
 BEH HILIC = mid pH (pH 5), 10 mM ammonium acetate with 0.02% acetic acid in 95:5 ACN:H<sub>2</sub>O  
 BEH Amide = high pH (pH 9), 10 mM ammonium acetate with 0.04% ammonium hydroxide 95:5 ACN:H<sub>2</sub>O  
 BEH Amide [Alternate] = mid pH (pH 5), 10 mM ammonium acetate with 0.02% acetic acid in 95:5 ACN:H<sub>2</sub>O

\*BEH represents either XBridge HPLC columns or ACQUITY UPLC BEH columns

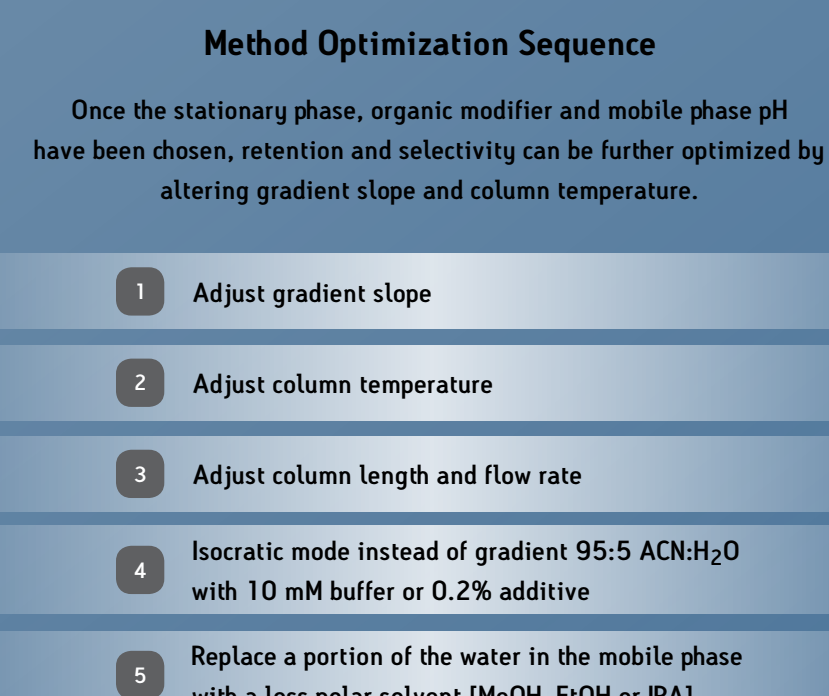

Evaluate result after each step. Stop after criteria for success has been met. Consider injection solvent (sample diluent) if poor peak shape/resolution.