


Conditions

Column: ACE 5 C18
Dimensions: 250 x 4.6 mm
Part Number: ACE-121-2546
Mobile Phase: MeCN/50 mM KH_2PO_4
pH 4.4 in H_2O (40:60 v/v)
Flow Rate: 1 mL/min
Temperature: Ambient
Detection: UV, 230 nm

Analytes

1. Phthalic acid
2. p-Hydroxybenzoic acid
3. Benzoic acid
4. Sorbic acid
5. Methyl paraben
6. Ethyl paraben
7. Propyl paraben

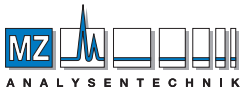

AUTHORIZED DISTRIBUTOR

MZ-Analysentechnik GmbH
Barcelona-Allee 17 • D-55129 Mainz
Tel +49 6131 880 96-0
Fax +49 6131 880 96-20
e-mail: info@mz-at.de
www.mz-at.de

