



## Conditions


Column:	ACE 5 HILIC-A
Dimensions:	150 x 4.6 mm
Part Number:	HILA-5-1546U
Mobile Phase:	12 mM ammonium formate pH 4.7 in MeCN/H <sub>2</sub> O (88:12 v/v)
Flow Rate:	1.5 mL/min
Injection:	2 μL
Temperature:	25 °C
Detection:	UV, 230 nm

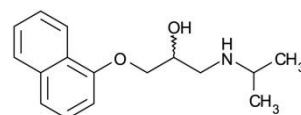


1. Propranolol

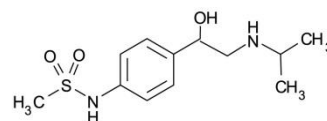


2. Sotalol

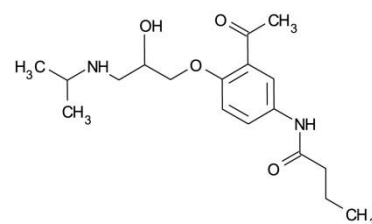

3. Acebutolol

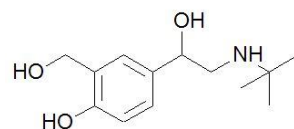

4. Salbutamol

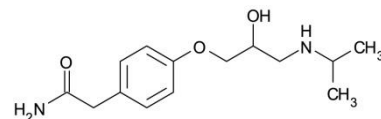

5. Atenolol

