



Analysis of Free Amino Acids

Data No. LA011-0000

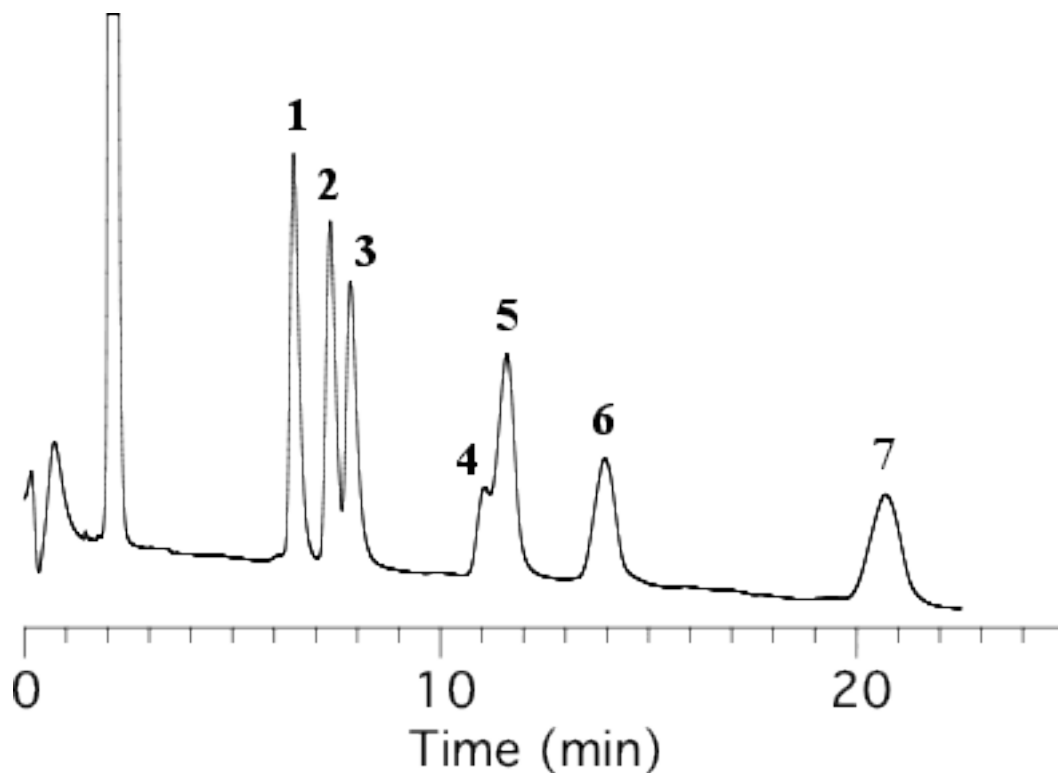

Column	: Inertsil Diol 5 µm 150 x 4.6 mm I.D.	Solutes	: 1. Phenylalanine 2. Methionine 3. Valine 4. DOPA (3,4- dihydroxyphenylalanine) 5. Alanine 6. Serine 7. Citrulline
Column Cat. No.:	: 5020-05645		
Eluent	: A) CH ₃ CN B) H ₂ O A/B = 85/15, w/w		
Flow rate	: 1.0 mL/min.		
Col. Temp.	: 40°C		
Detection	: RI		
Injection Vol.	:		
Sample	:		

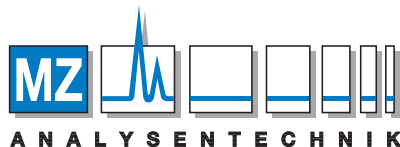

AUTHORIZED DISTRIBUTOR

MZ-Analysentechnik GmbH, Barcelona-Allee 17 • D-55129 Mainz
Tel +49 6131 880 96-0, Fax +49 6131 880 96-20
e-mail: info@mz-at.de, www.mz-at.de