


## Analysis of water-soluble vitamins

Data No. LA507-0000

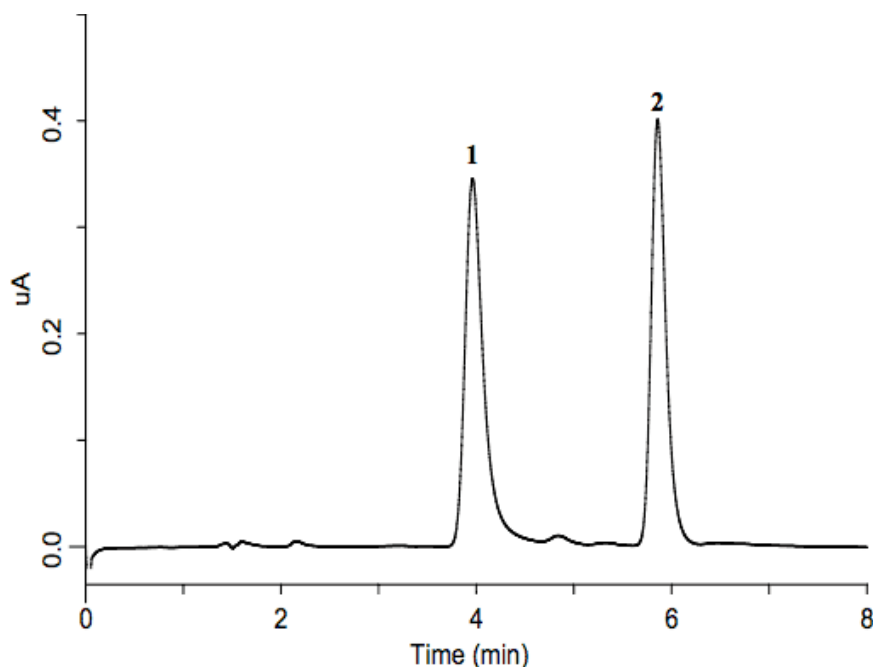

System : GL-7400 system

Column : Inertsil ODS-SP (5  $\mu$ m, 150  $\times$  4.6 mm I.D.)

Column Cat. No. :

Eluent : A) CH<sub>3</sub>CN

B) 0.1 % H<sub>3</sub>PO<sub>4</sub>

A/B = 18/82, v/v

Flow rate : 1.0 mL/min

Col. Temp. : 40 °C


Detection : ECD 1600mV vs. Ag/AgCl (ED 703 ECD Detector, Diamond)

Injection Vol. : 10  $\mu$ L

Sample : water-soluble vitamins

Solutes : 1. Cyanocobalamin (2  $\mu$ g/mL)

2. Biotin (2  $\mu$ g/mL)

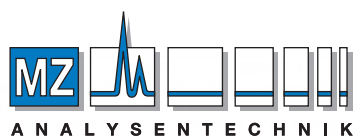

#### AUTHORIZED DISTRIBUTOR

MZ-Analysentechnik GmbH, Barcelona-Allee 17• D-55129 Mainz

Tel +49 6131 880 96-0, Fax +49 6131 880 96-20

e-mail: info@mz-at.de, www.mz-at.de