


Application Number 107920
 Method HPLC
 Results

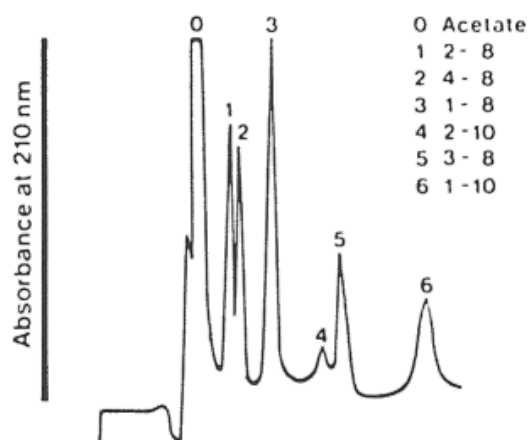

Legend:

Baseline resolution of 100 ng standard ANG-(1-8) octapeptide (peak 3) from other ANG peptides by isocratic HPLC

Substances angiotensin (ANG) peptides

Product(s)

Phase	REF	Webshop
NUCLEOSIL C18	720023.46	Shop now

Matrix plasma

Sample(s)

Conditions

Eluent: methanol - 0.085 % phosphoric acid (33.5:66.5, v/v)
 Flow rate: 1 mL/min
 Temperature: 45 °C

Detection

UV, 210 nm

Note

Author

Nussberger, J.


Source


1984 Blood Pressure Council Suppl. I Hypertension
7 (1985) No. 3, 1 - 7

Disclaimer

Concerning questions to literature applications please contact the cited author/s. Chemicals mentioned in applications are not manufactured or sold by MACHEREY-NAGEL. Applications mentioned in MACHEREY-NAGEL literature or databases are provided for informational purposes only. MACHEREY-NAGEL does not warrant that all applications have been tested in MACHEREY-NAGEL laboratories using MACHEREY-NAGEL products. MACHEREY-NAGEL does not warrant the correctness of the applications.

Share via Mail

 Print View

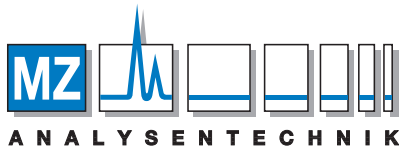

AUTHORIZED DISTRIBUTOR

MZ-Analysentechnik GmbH, Barcelona-Allee 17 • D-55129 Mainz

Tel +49 6131 880 96-0, Fax +49 6131 880 96-20

e-mail: info@mz-at.de, www.mz-at.de