

1 Glycolic acid

2 Formic acid

3 Acetic acid

HPLC Conditions


Column CAPCELL PAK C18 MG S5 : 2.0 mm i.d.×250 mm

Mobile phase 0.1vol% H₃PO₄, H₂O / CH₃CN = 97.5 / 2.5 (v/v)

Flow rate 200 µL/min

Temperature 35 °C

Detect UV 210 nm

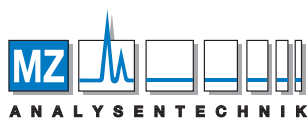

AUTHORIZED DISTRIBUTOR

MZ-Analystechnik GmbH, Barcelona-Allee 17 • D-55129 Mainz

Tel +49 6131 880 96-0, Fax +49 6131 880 96-20

e-mail: info@mz-at.de, www.mz-at.de