

- | | | |
|----------------------|------------------|------------------|
| 1 PTH-Aspartic acid | 2 PTH-Asparagine | 3 PTH-Threonine |
| 4 PTH-Glutamic acid | 5 PTH-Glycine | 6 PTH-Alanine |
| 7 PTH-Tyrosine | 8 PTH-Methionine | 9 PTH-Tryptophan |
| 10 PTH-Phenylalanine | 11 PTH-Leucine | |

HPLC Conditions


Column CAPCELL PAK C18 UG120 S5 : 1.5 mm i.d.×250 mm

Mobile phase (A) 20mmol/L CH₃COONa (pH4.0)
 (B) CH₃CN
 B 15% (0.0min) - 20% (9.0min) - 45% (17.0min) - 45% (30.0min) Gradient

Flow rate 100 µL/min

Temperature 35 °C

Detect UV 254 nm

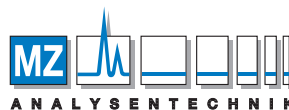

AUTHORIZED DISTRIBUTOR

MZ-Analysentechnik GmbH
 Barcelona-Allee 17 • D-55129 Mainz
 Tel +49 6131 880 96-0
 Fax +49 6131 880 96-20
 e-mail: info@mz-at.de
 www.mz-at.de