

Data code:


T01016

Application Title:

Separation of acrylamide by HILIC

Compound:

acrylamide

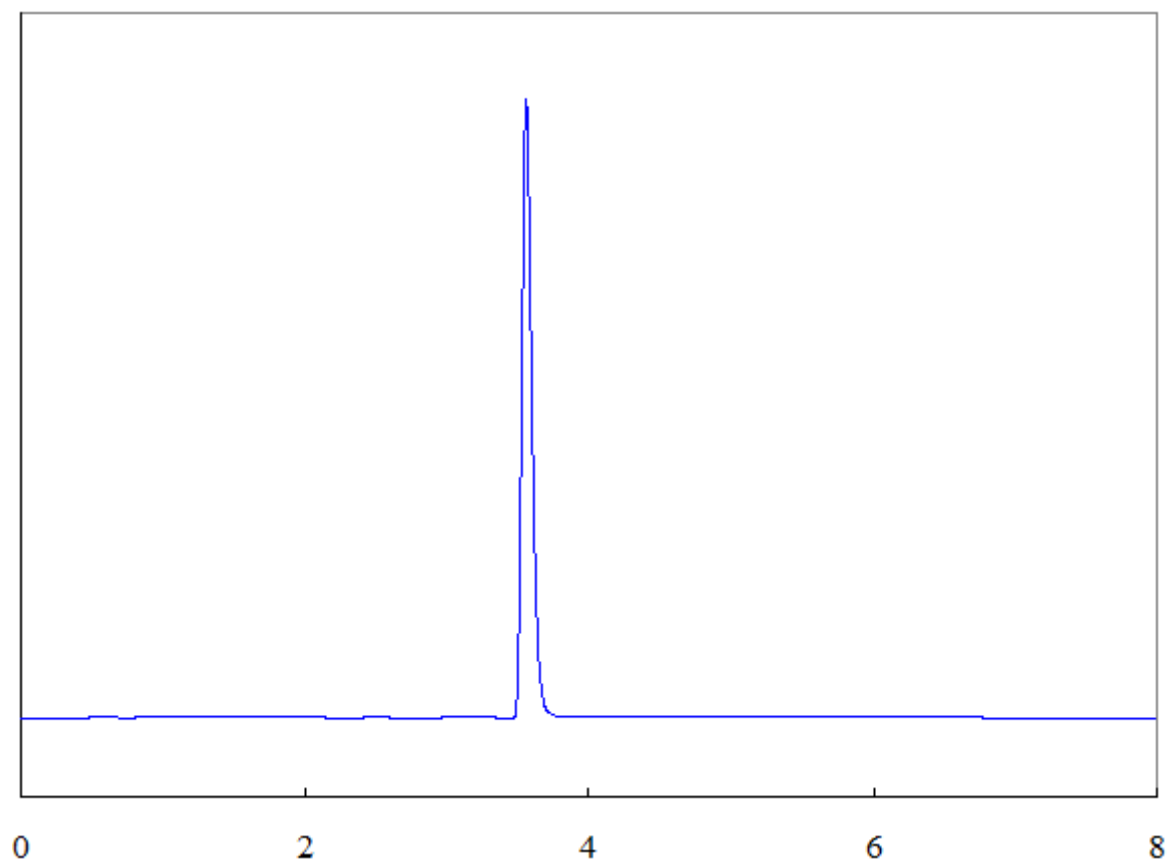

Chromatographic conditions

Column: TSKgelAmide-80 5 μ m 4.6mm \times 25cm, P/N 13071, 4.6mm \times 25cm

Mobile Phase: H₂O/CH₃CN = 20/80

Flow Rate: 1.0mL/min

Detection: UV/VIS 254nm

Temperature: 40° C

Injection Volume: 20 μ L

Title: Separation of acrylamide by HILIC